

Turbidity

An interesting parameter

Brief Overview, Measurement Tips & Tricks



Turbidity – Driven by Application

- **Drinking Water Applications**

- Indicator for water purity
- Quality Control of ponds, reservoirs and water works
- Control of (over) use
- Drinking water in Food & Beverage industries
- Reference measurements for process instrumentation

- **Waste Water Applications**

Reference measurements for process instrumentation

- **Industrial Process Control, Production, QC of incoming goods**

- Filter loads: avoiding breakage of filter
- QC of fuel quality: income control in motor development, ...
- Food & Beverage: production like juice, washing water of process...
- Commercial Fish-farming
- Cell culture growth, fermenter
-

- **Environmental monitoring: Surface Water Bodies**



Measuring Principles

Standard methods & application driven




a xylem brand

Principles of Turbidity Measurement


Turbidity is an optical impression!

Thus, turbidity is not reflecting a clearly defined measure like e.g. electrochemical or photometric parameters, but is an effect caused by scattering and absorbance of light instead:

- Solid particles cause light scattering, depending on size, shape and colour
- Angles of light scattering are varying with size and particle count: forward, sideward and backscattering.
- Attenuation of light is increasing by multiple scattering and absorbance with growing particle count (= higher turbidity) => Scattering intensity is non-linear within the complete range

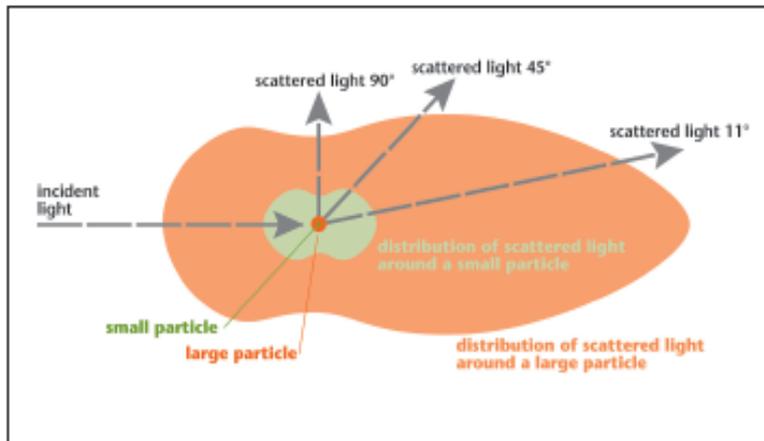

Scattering effects are dependent on particle size: green is a small, orange is a big particle

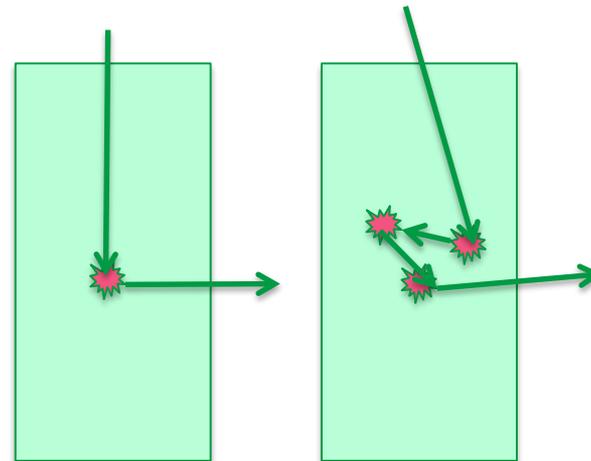


The more particles the more multiple scattering between particles occur

Principles of Turbidity Measurement

What do scattering effects mean for readings and reproducibility?

Imagine a piece of paper in the light: being exposed full formate or from the side

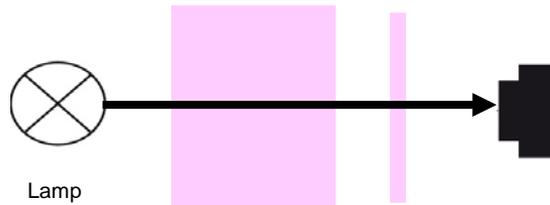

- Particles are not stable and therefore are changing their orientation in the solution
- With changing orientation of an unevenly shaped particle in the light, the scattering effect will be varying:

This results for the measurement in:


- „Oscillating“ readings around an average value
- Higher tolerances / variability of the system, measurement range dependent
- Resolution decreases with higher measurement range due to interactions
- In comparison to other measurement procedures, a real sample does not show a „stable“ measurement system during measurement.

The three main procedures

Nephelometric - Scattering at 90° angle

Best suitable measurement for drinking water applications
Standard method of DIN EN ISO and US EPA


- Light beam passing through a cell
- Scattering detection with a 90 detector orientation
- Light passing through the cell is being ignored

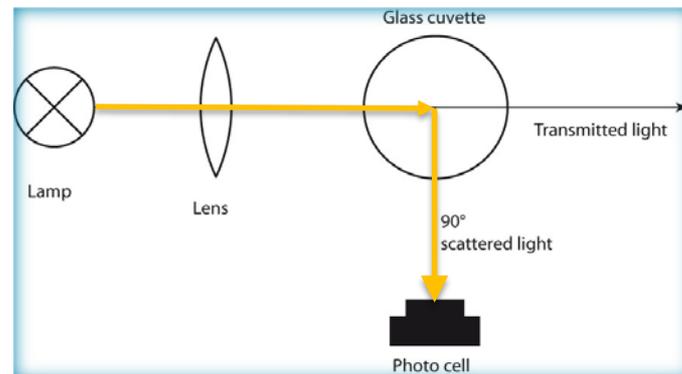

Transmission measurement at 180°

- Transmisson is defined by light intensity after/before vial
- The detector is placed at 180°
- Photometric measurement is performed via transmission:


This is the principle of photometric turbidity correction for some parameters like COD or being used for photometrical QC of e.g. liquids like wall paint.

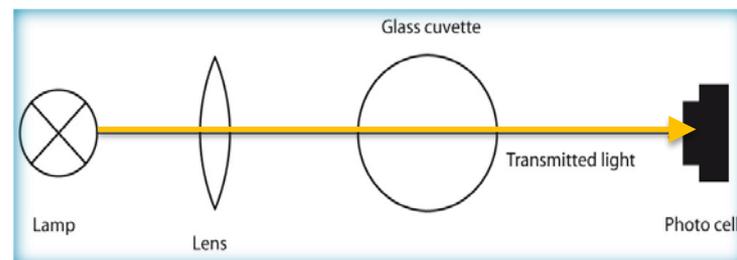

Three main procedures

Ratio measurement – combining various detection angles

- Important for solutions with high turbidity values >1500- 2000 NTU to optimize detection of multiple (back) scattering effects between particles.
- Measuring different (scattering) angles in ratio to transmission being evaluated algorithmically
The more particles are dissolved, the more interacting effects are given. This may lead to an even lower measurement value with the 90° detection depending on particle distribution. With the ratio procedure a detection up to 10.000 NTU is possible (taking a high tolerance into account).
- **No standardized ratio procedure => no comparable results between different instruments**
Application specific standards, e.g. brewery industry with defined application optics (90°, 11°)

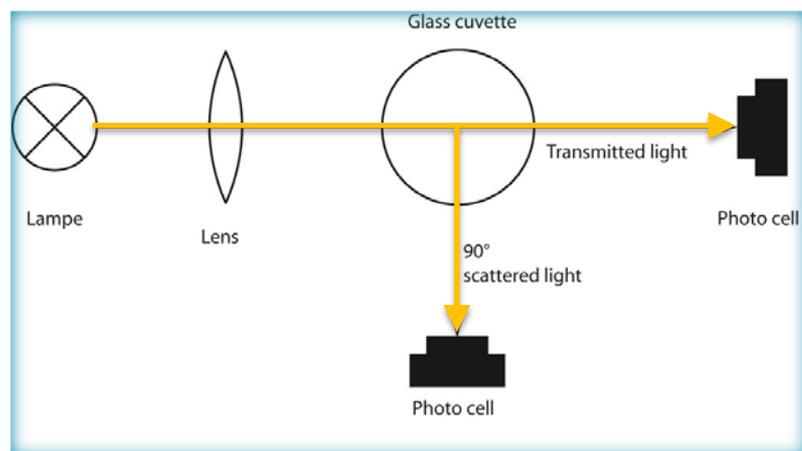

Attention in water treatment referencing process with lab meters:

Results from lab instruments with ratio from 40 NTU onwards are not comparable to nephelometric results above 40 NTU! Ratio is a non-standardized but manufacturer dependent optical system. Values should have any „Ratio“ mark to be identified as such.

IR or Tungsten Halogene Lamp?

Different Light Sources

IR at 860 nm

- fulfills **DIN EN ISO 7027-1** requirements => European standard for drinking water
- tech. requirements of optical system strictly defined:
 - 860 nm \pm 30 nm
 - - Detector at $90^\circ \pm 2.5^\circ$, aperture angle 20-30°
- IR eliminates colour effects of solutions, less influence of ambient straylight

White light - Tungsten Halogene

- fulfills **US EPA 180.1** requirements => US regulations for drinking water:
- Technical specification is much wider
 - light source between 400-600 nm, filament temperature 2200 – 3000 K
 - Detector at $90^\circ \pm 30^\circ$, aperture $< 30^\circ$
- In the given wavelength range, color influences by absorbance are mostly disturbing
- Due to smaller wavelength, smaller particles can be detected

=> Suitable optical systems are depending on standard methods, industries and applications

What is the best choice?

If not defined by requirements of standard methods, the application is important for the best choice:

Particle size and count + coloration of samples!

Measuring task \ Optics	IR 90° 860 nm	T 90° 400-600 nm	Forward scattering (IR or T) (e.g. 11-45°)	Ratio (IR or T)
Low Turbidity Level e.g. Drinking Water <1 NTU	x	x		
High Turbidity Level > 1000-2000 NTU				X (cannot be compared to nephelometric results)
Colored samples	x			x
Small Particles		x		
Big Particles, e.g. cells	x		x	


Influence of color – Example Fruit Juices

Samples of fruit juices, turb and clear	Comparison IR vs T			
	Turb 430 IR	Turb 750 IR	Turb 430 T	Turb 750 T
Beetroot	220	222	56	50
Cabbage (turb)	789	802	695	683
Carot	556	556	112	115
Vegetable	764	770	138	128
Blueberry	708	746	31,2	25,5
Banana (turb)	996	983	350	354
Cherry-Banana	OFL	OFL	253	236
Veg/Tomato	583	608	85	77
Apple clear	0,52	0,43	0,08	<0,01

IR light at 830-860 nm is out of absorbance range of colored solutions:

- ⇒ The measurement results are mainly defined from turbidity alone: IFU is recommending IR for measurement in juice.
- ⇒ IR optics allow less variation in lamp and layout: Measurement results will show mostly variability between different instruments.

Tungsten lamp with a variety of 400-600 nm wavelength AND color temperature is in the absorbance range of colored solution:

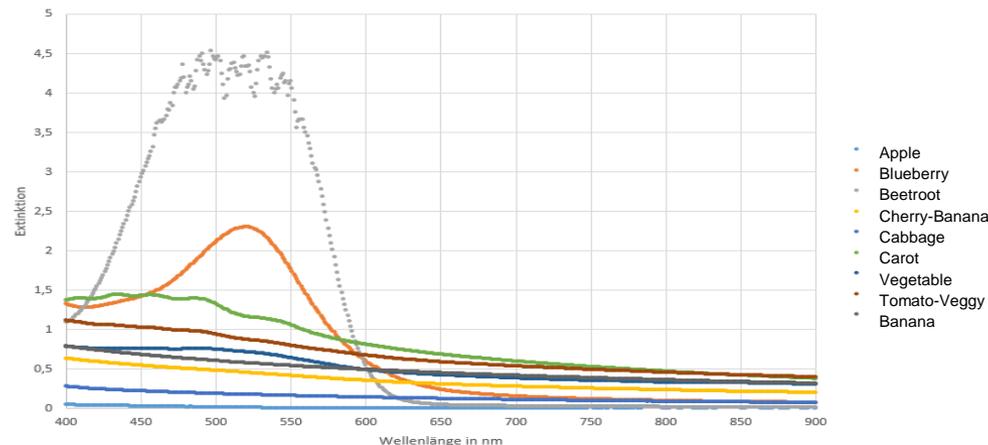

- ⇒ The spectral absorbance in the visible range (see spectra) influences the turbidity measurement extremely.
- ⇒ The designs (lamp, detector) are different: coloration and color intensity increase the variability in results between different instrument models (from different brands).
- ⇒ With high turbidity the light – other than with clear solutions - is probably not pathing through the complete vial path length.

Photometric spectra with photoLab® 7100 VIS, 1:20 dilution

More factors are color intensity of sample, lamp color temperature of lamp and detectors (not investigated here)

Calibration & Standards

Turbidimetry has natural limitations



a xylem brand

Calibration

Brief History of Calibration:


The Jackson Candle Turbidimeter explains the calibration idea of turbidity meas

Attempting a “real world sample image” with differently shaped particles

- First standards from diatoms in earth (= **Kieselguhr**)


Diatoms are algae with silica housing and **multiple** shape and size
Prepared from a 1000 ppm stock solution = suspended silica dilution steps

*) Source Image: <https://de.wikipedia.org/wiki/Kieselguhr#/media/Datei:Diatoms.png>
Images courtesy of Mary Ann Tiffany, San Diego State University. - Bradbury J:
Nature's Nanotechnologists: Unveiling the Secrets of Diatoms.
PLoS Biol 2/10/2004: e306. [doi:10.1371/journal.pbio.0020306](https://doi.org/10.1371/journal.pbio.0020306)

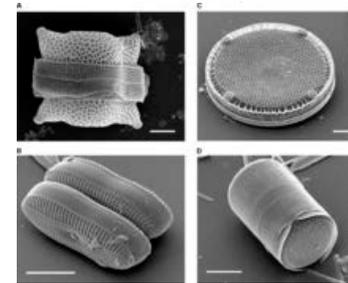

Various diatoms,
electronmicroscope*)

- Poured from top view in glass tube with candle below
- With increasing particles light turns gradually to a uniform shimmer and final opaque solution. These grades have been marked
- Unit = JTU, Jackson Turbidity Unit, limitation of this „transmission“ < 25 JTU

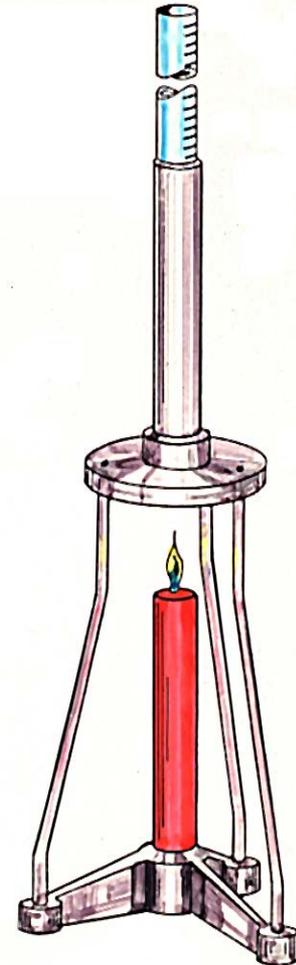


Calibration Yesterday and Today

Formazine in the 1950ties: **best image of real samples**

Primary standards DIN EN ISO and US EPA, characterized by a high tolerance need


- Hazardous, cancerogene (raw) material
- Stock solution of 4000 NTU, variation of 5-10% (manufacturer depende depending on raw material and manufacturer, additional dilution error adds
- Freshly prepared dilutions for use within 24 hrs only, 2 hrs for small values!
- Particle instable in shape and size by clogging and deterioration of particles
=> **inhomogenous and varying dispersion (even for stabilized formazine)**

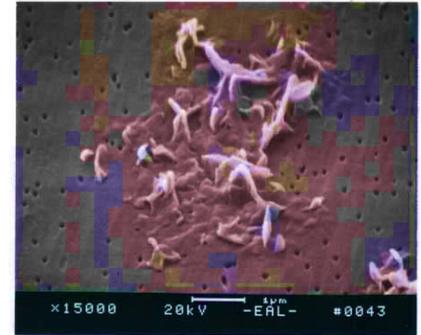

AMCO Clear® Standards and other polymers: **Offering best and stable real sample image today**

Polymer microspheres in ultrapure water are accepted to be primary standard acc. US EPA, secondary standard acc. to DIN EN ISO, highly precise and longterm stable


- Not hazardous
- **1% Lot-to-Lot accuracy** in production for size, shape
=> stable and homogenous distribution in particle composition
- Traceable to Formazine
- Designed and optimized for each optical system:
=> **Highest precision** reached with specific IR and T calibration standards
=> Only freshly prepared Formazine standards from **one and the same standard for an instrument comparison study**

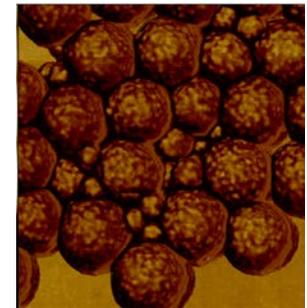


Calibration: Comparing Standards

AMCO Clear® Standards (1) vs Formazine (2)

AMCO Clear® Standards	Formazine / stabilized Formazine
Non hazardous	Formazine raw material is rated carcinogenic and hazardous
Defined polymer beads, unchanged composition, no dilution error 1% Lot to Lot accuracy	Formazine stock solution plus dilution steps lead to 10- 20 % deviation between different manufacturers
Stable particle distribution and particle size, can be optimized for each optical system	Drift of particle size and shape => values drift accordingly. Stabilized Formazine standard need "Re-Suspension" => Change in size and distribution
Easy to dispose	Formazine to be disposed separtely
Longterm stability min. 15 months in bottles, min. 3 months in clean vials	Formazine Standards change within 24 hours, stabilized solution also drift (see above)

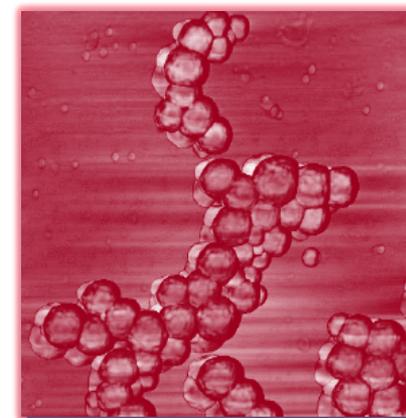

Polymers

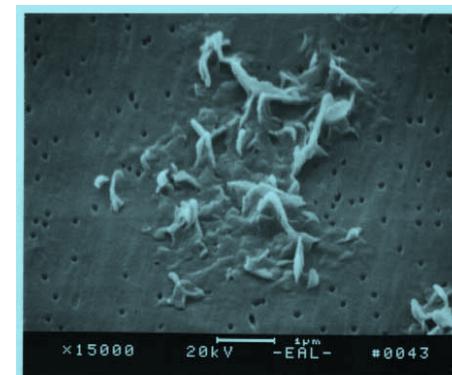

Formazine

How to achieve best results?



a xylem brand

Turb[®] 750/430: Simply „all right“ - with „IRPC“

„Stable“ is, what is correct!

Precise results in turbidity measurement mean:

Measuring moving, instable particles reliably with **IRPC**

Intelligent **R**eproducibility- and **P**lausibility **C**heck

=> Internal multiple readings per measurement

=> Evaluation of these values

⇒ Elimination of outliers

⇒ **Stable data output of the evaluated result as a stable and reproducible reading!**

Turbidity	
[IRPC]	[AQA]
157.0	FNU NTU
[3-P StdCAL]	
15.01.20 15:12	

Take limits of the measurable into account

⇒ Turbidimetry with formazine showing 10% tolerances and more has to be realistic in terms of expectation: depending on measurement range 2-3% of the value at higher ranges is **very good!**

⇒ Higher tolerances are no problem: turbidity values are giving a value corridor of „being OK“ for quality and contaminant rating

⇒ Resolution of 0,02 NTU in low range can only be achieved with best practice and correct handling! Smaller resolution 0.00x are simply digital. Water scatters with approx 0,02 NTU/FNU itself.

Best measuring procedure: preparation

Excluding error sources

- Dirt and Dust: Clean vials before use, no fuzz in calibration standards, clean vial from finger prints with optical
- Glass inhomogeneities for better reproducibility: Indexing (see below)
- Scratches excluded by indexing for suitable position of vial (see below)

Prepare vials by indexing for longterm-use

Recommendation in regulations (e.g. US EPA)

- Place vial, press START and slightly turn to the right and left
- Release START when a „stable“ value range: meaning no big jumps up or down, which would indicate a scratch*)
- Mark position on cap label

*) Tolerance for value deviations:

± 0,02 NTU below 1 NTU

± 2-3 % for higher values

Then measuring routine is fast, precise, comfortable

- Place prepared = marked vial
- Press and release START
- Value will be stored, or store with ID



Best measuring procedure

Advantage of Indexing versus pre-marked vials

is the optimized path of light through the vial without impacts

- Scratches and glass inhomogenities are always excluded- in contrary to pre-marked vials - with indexing procedure and allow longterm use
- Big value jumps during indexing indicate scratches or glass faults => indexing allows different positions to lows value for scratch-free light path
- Scratches in light path result in wrong turbidity values, mostly higher => pre-marked vials result eventually in less accurate values
- **Silicon oil** to mask scratches eventually has a differing index and **attracts dust and dirt** on vials => misleading results
- **For single use of vials:**
Search best position by aligning the vial with lowest value during measurement procedure as displayed.



How to achieve optimal results in drinking water?

Assuring purity and cleanliness for human consumption!

- ⇒ Value ≤ 1 NTU, since particles offer ground for bacterial/viral growths
- ⇒ This value is very low considering the limits of this measurement procedure
- ⇒ Accuracy of result and assurance require a maximum accuracy during calibration and measurement procedure
- ⇒ Dirt and Dust: Clean vials before use, no fuzz entry into calibration solutions, wipe off fingerprints with optical cloth.
- ⇒ Avoid glass inhomogeneities and scratches by indexing



Calibration of the 0,02 NTU/FNU standard

- Insert vials, slightly turn to the left and right
- Release START when value is lowest respectively in a „stable range“ w/o big jumps of value, which might indicate a scratch (see manual)
- Mark position on cap label

Use same vial for measurement of samples and 0.02 NTU/FNU standard

Sample measurement in the same vial as 0.02 NTU standard before empty vial, rinse with sample, fill (next) sample, measure at aligned marking position each time

The 0.02 NTU standard (with preservative against bacterial contamination) can be substituted by fresh, clean and particle-free de-ionized water.



GLP – The daily proof

A classical chicken-and-egg issue: What is drifting - standard or instrument?

- In laboratories meters drift less than in field situations due to stable ambient condition.
- But: for good laboratory practice, a proof standard before measuring can be used.

Tips for practise of a proof standards for low values

- ✓ Prepare a dilution out of a 10 NTU Polymer-Standards to a nominal value of 0,3x to 0,6x NTU
- ✓ Daily measurement and documentation of the value in a tolerance range of $\pm 10\%$
- ✓ Does the proof exceed $\pm 10\%$: Perform a calibration and prepare a new proof standard

Date	xx.xx.xx	xx.xx.xx									
												> 10%
0,36												
0,35						X						
0,34		X							X			
0,33±0,033	X		X		X		X	X		X		
0,32				X							X	
0,31												
0,30												
												X

Avoiding Impacts

- **Correct sampling and fast measurement:**
 - Bring sample to room temperature, homogenisation
 - Do not shake vigorously (bubbles => wait)
 - Re-suspension by gentle movement
 - Fast settling samples instantly, average values
 - No condensation on vial
- **Impacts:**
 - Bubbles in the sampl
 - Temperature change
 - Condensation (Temp!)
 - Finger prints!
 - Scratches => Indexing
 - Dirt and fuzz in vials
- **The correct optical choice to match the application tasks**



Summary

- Turbidity is an indicator parameter: E.g in drinking water, particles stand for a platform of bacterial growth, meaning a potential hazardous bacterial contamination.
- Turbidity measurement principles are not comparable to physical or electro-/chemical measurement principles, where clearly defined quantities or concentrations of substances are reacting and can be analyzed subsequently. E.g. pH, Ammonium, Chlorine, COD...
- For turbidity measurement, various optical systems can be selected matching different application requirements and purposes: E.g. does color rarely influence the measurement with IR lamp negatively in comparison to tungsten lamp.
- Formazine – beside using harmful raw material – is less accurate than modern polymer particles with a production accuracy of +/- 1% like the well proven AMCO Clear® standards.
- The comparison of measurement results from different instrument models is only given when all instruments are calibrated with one and the same freshly prepared formazine standards.
- For drinking water and low values ≤ 1 NTU the accuracy can be increased by following an improved calibration and measurement handling together with instrument specific or other polymers.

For more information:

www.xylemanalytics.com

Technical Support

Email: TechInfo.WTW@xylem.com

Phone: ++49 881 183-321



a xylem brand